

MR Brain Routine IVcon wo/w

MR Brain Routine
IVcon wo/w


IMG15436

Acquisition Information


	Plane	Sequence Name	Thickness - Spacing	FOV
1	Axial	DWI (B=0,1000)	4 mm - 1 mm	220-260 mm
2	Axial	T2 TSE	4 mm - 1 mm	Approx. 200 mm
3	Sagittal	T1 FLASH	4 mm - 1 mm	Approx. 200 mm
4	Axial	T1 TSE	4 mm - 1 mm	Approx. 200 mm
Power Inject Gadavist @ >1.5 mL/sec				
5	Axial	T2 FS FLAIR TSE	4 mm - 1 mm	Approx. 200 mm
6	Axial	T1 TSE	4 mm - 1 mm	Approx. 200 mm
7	Axial	3D SWI	4 mm - 1 mm	Approx. 200 mm

	Plane	Sequence Name	Thickness - Spacing	FOV
8	Coronal	3D T1 MPRAGE Iso	1 mm - 0 mm	Approx. 200 mm
9	Coronal	T2 TSE	4 mm - 1 mm	Approx. 200 mm

Reference Images

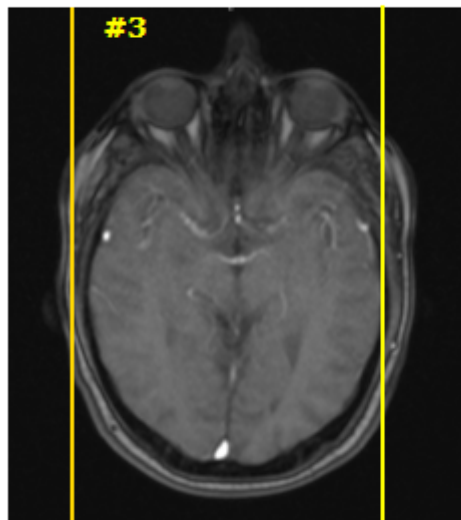

Align series #1,2,4-7 parallel with subcallosal plane

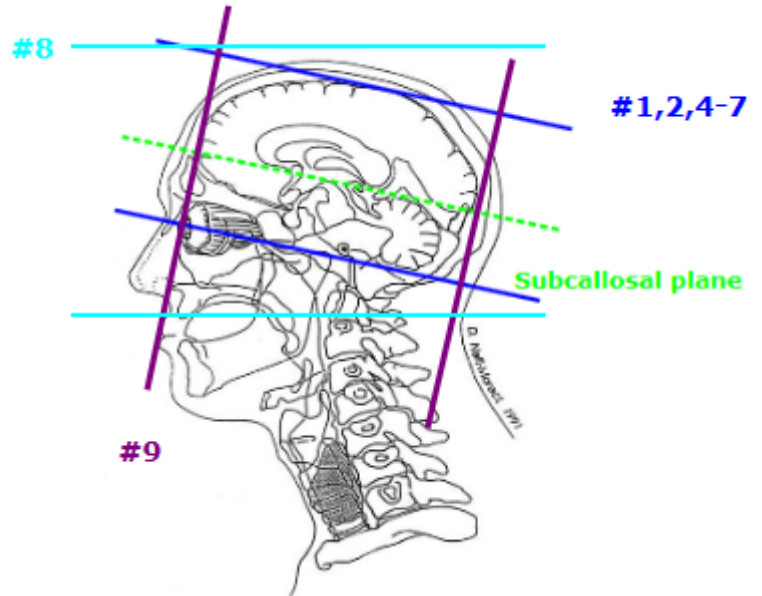

Special Instructions

*Ensure phase encoding is A>P

** Subtract Pre and Post; FOV, Slice Position, and number of slices must match for subtraction

*** must be at least 2 concatenations/covers/acquisitions

Indications

Revision #2
Created 2026-06-23 14:41:20 UTC by Justin Brown
Updated 2026-06-23 15:22:05 UTC by Justin Brown